

## The key to accurate positioning of Dr Suni Sensor Holders.

XCP-DS ® System for Suni Digital Sensors make it possible to position digital sensors perfectly every time.



- Securely holds sensor in place while helping to precisely align the x-ray cone
- Sensor holders (baskets) assure perfect coordination with Dr. Suni Sensors
- Standardized aiming rings, arms and bite-pieces simplify reordering
- Ergonomically designed for greatest patient comfort
- Fully color-coded for easy assembly
  - Blue = Anterior
  - Yellow = Posterior
  - Red = Bitewing
  - Green = Endo
- Available in comprehensive and specialized kits



XCP-DS® Replacement Parts

The following components are common to all XCP-DS systems.


Sold individually

- 54-0865 **Anterior** Aiming Ring
- 54-0860 **Posterior** Aiming Ring
- 54-0934 **Bitewing** Aiming Ring
- 55-0598 **Endo** Aiming Ring



Sold individually

- 54-0857 **Anterior** Arm
- 54-0858 **Posterior** Arm
- 54-0927 **Bitewing** Arm
- 55-0597 **Endo** Arm

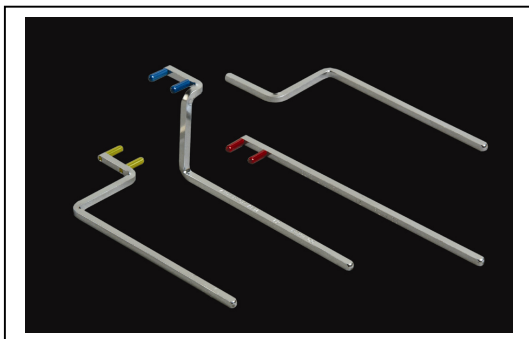




3 Bite pieces per package

55-0567 **Anterior** Bite Piece

55-0568 **Posterior** Bite Piece

55-0570 **Bitewing** Bite Piece

55-0562 **Endo** Bite Piece

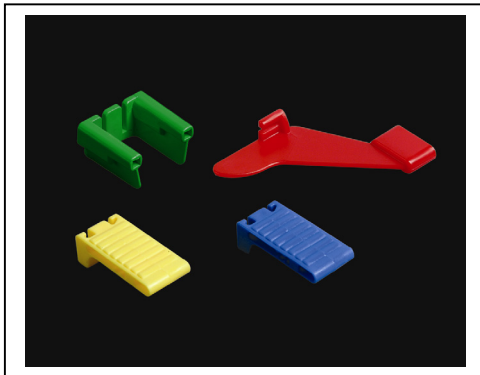

#### **55-1604 Combo Kit**

Kit includes:

4 Arms: 1 Anterior, 2 Posterior, 1 Bitewing

8 Bite Pieces: 2 Anterior, 4 Posterior, 2 Bitewing

4 Aiming Rings: 1 Anterior, 2 Posterior, 1 Bitewing

20 Baskets (Holders): 5 #2 Horizontal, 5 #2 Vertical, 5 #0 Horizontal, 5 #1 Vertical

Instructional Poster

#### **55-1600 #2 Kit Complete**

Kit includes:

4 Arms: 1 Anterior, 2 Posterior, 1 Bitewing

8 Bite Pieces: 2 Anterior, 4 Posterior, 2 Bitewing

4 Aiming Rings: 1 Anterior, 2 Posterior, 1 Bitewing

10 Baskets (Holders): 5 #2 Horizontal, 5 #2 Vertical

Instructional Poster

#### **55-1601 #2 Kit**

Kit includes:

8 Bite Pieces: 2 Anterior, 4 Posterior, 2 Bitewing

10 Baskets (Holders): 5 #2 Horizontal, 5 #2 Vertical

Instructional Poster

**55-1602 #0/1 Kit Complete**

Kit includes:

4 Arms: 1 Anterior, 2 Posterior, 1 Bitewing

8 Bite Pieces: 2 Anterior, 4 Posterior, 2 Bitewing

4 Aiming Rings: 1 Anterior, 2 Posterior, 1 Bitewing

10 Baskets (Holders): 5 #0 Horizontal, 5 #1 Vertical

Instructional Poster

**55-1603 #0/1 Kit**

Kit includes:

8 Bite Pieces: 2 Anterior, 4 Posterior, 2 Bitewing

10 Baskets: 5 #0 Horizontal, 5 #1 Vertical

Instructional Poster


**Replacement Parts**

**55-0589** #2H Basket, 3 pcs.

**55-0590** #2V Basket, 3 pcs.

**55-0591** #0V Basket, 3 pcs.

**55-0592** #1V Basket, 3 pcs.

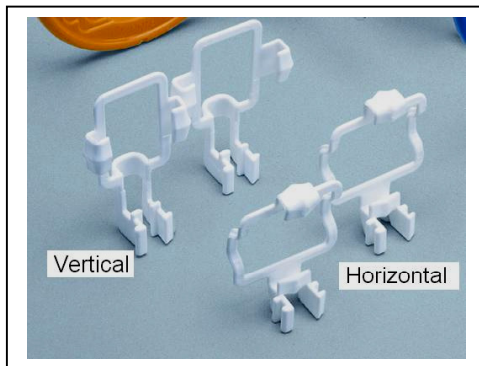

**5554** Instructional Poster

XCP-DS® Dr. Suni Sensor Holders

**Placing the sensor in the holder**

1. Pull a silicone sensor sleeve over the sensor and stretch it out along the cable.
2. Slip one edge of the sensor into the holder and gently pull the cable through the base. Be sure the flat surface of the sensor faces the gripping edges and the back surface sits in the 'window' of the holder. (A)
3. Snap the second edge of the holder over the sensor. (B)
4. Slide the holder into the bite piece. (C)
5. Attach arm and ring.

